



## CORPECTOMY

NEXXT MATRIXX® 3D Printed Porous Titanium



With INERTIA® Pedicle Screw System

With STRUXXURE® Anterior Cervical Plate

**NEXXT**  
Spine

317-436-7801

INFO@NEXXTSPINE.COM

14425 BERGEN BLVD, STE B

NOBLESVILLE, IN 46060


70-069E, Rev A



## NEXXT MATRIXX®


3D Printed porous titanium

10X

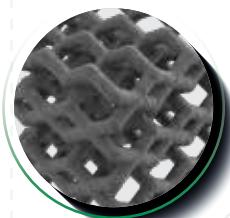

Systematic Titanium  
PORES

300X

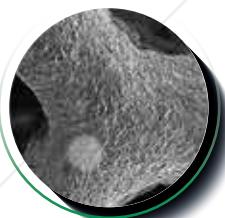

Uncompromising  
MACROSURFACE


10,000X



7 $\mu$  Surface  
MICROSURFACE

## CORPECTOMY

Total implant lattice integration

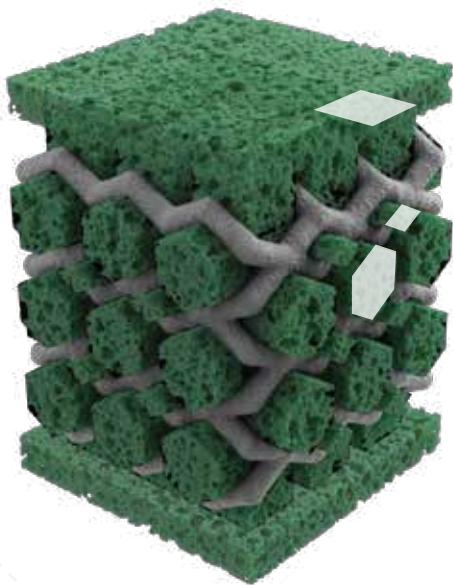

*Image represents potential bony ingrowth*

### TI PORES

- NEXXT MATRIXX® exhibits three pore sizes of 300, 500, and 700 $\mu$ .
- Minimized titanium material resulting in a 75% open porous architecture.


### SURFACE

- Nexxt Spine has developed a proprietary, residue free, micro-roughening process that creates a highly cohesive 7 $\mu$  roughened topography.
- Due to the micro-roughened porous structure of the NEXXT MATRIXX® titanium, the implants exhibit up to 4X more surface area for bone apposition and potential bony integration than conventional spinal implants.




## FEATURES


Enabling new corpectomy implant possibilities

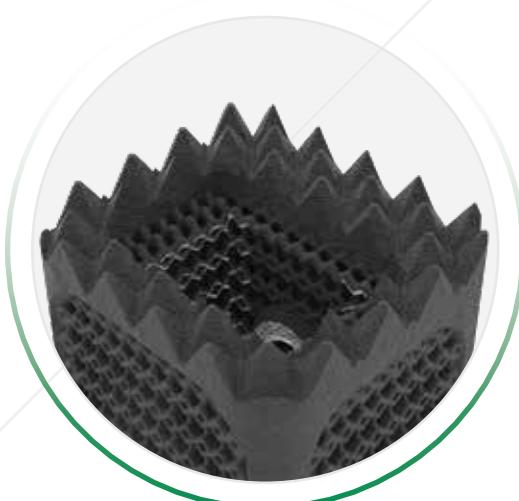

Aggressive teeth for secure placement

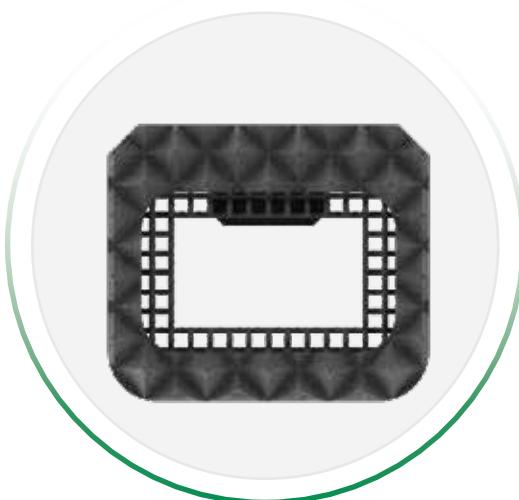


Entire implant graft support

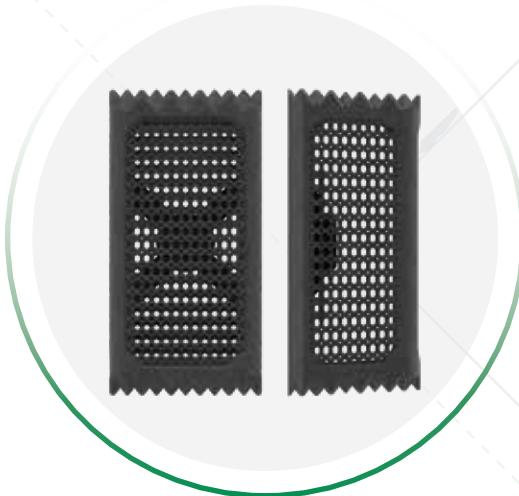

Complete A/P + lateral visualization

## INSTRUMENTATION

For intraoperative optimization

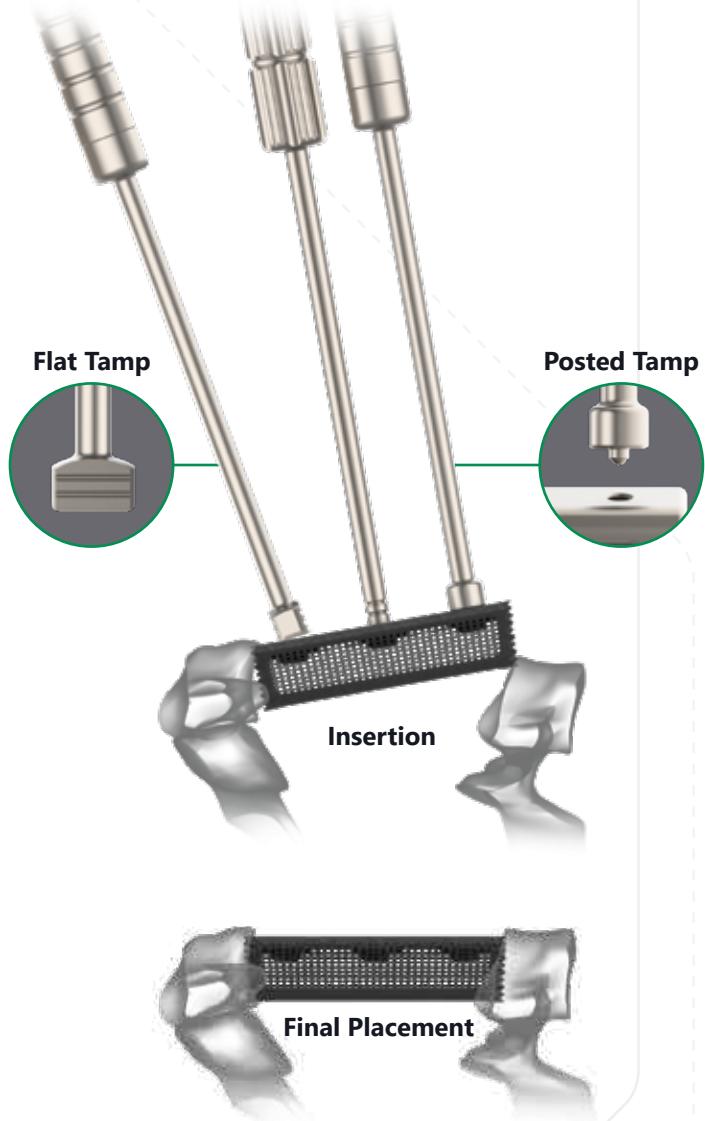

## IMPLANT OFFERING

Accommodation of anatomy

### 12D x 14Wmm 0°

12Dx14Wx16-40Hmm 0° 2mm increments

### 12D x 14Wmm 6°

12Dx14Wx11-30Hmm 6° 1mm increments

12Dx14Wx32-50Hmm 6° 2mm increments

### 14D x 16Wmm 0°

14Dx16Wx16-40Hmm 0° 2mm increments

### 14D x 16Wmm 6°

14Dx16Wx11-30Hmm 6° 1mm increments

14Dx16Wx32-50Hmm 6° 2mm increments