


**NEXXT MATRIXX®** Technology



One-Step Turn Lock

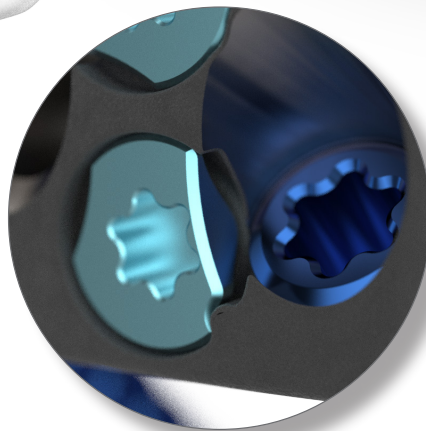

Self-Guided Instruments




Nexxt Spine, LLC  
14425 Bergen Blvd, Suite B  
Noblesville, IN 46060  
(317) 436-7801  
[Info@NexxtSpine.com](mailto:Info@NexxtSpine.com)

70-080E, Rev C




10X

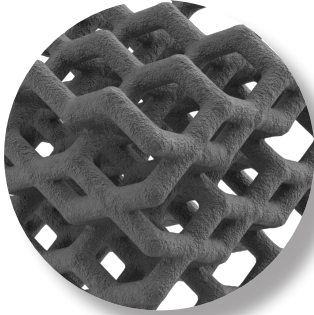

Systematic Titanium  
**PORES**

300X

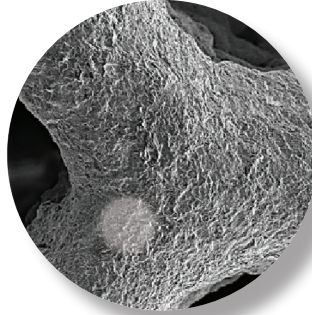

Uncompromising  
**MACROSURFACE**


10,000X



7µm Surface  
**MICROSURFACE**

### Pillars of NEXXT MATRIXX® Technology:

1. Varied pore array of 300, 500, and 700µm designed to support vascularization and osteogenesis.<sup>1,4,5</sup>
2. 7µm surface roughness designed to increase osteoblast differentiation, production of angiogenic factors, and surface osteointegration.<sup>2,3,6</sup>
3. 75% porous, open titanium architecture developed for greater surface area and nutrient exchange, leading to increased volume for potential bony in-growth.<sup>4,5,6</sup>
4. Modulus of elasticity engineered to be comparable to PEEK devices leading to a more physiological product.<sup>6</sup>
5. 700µm A/P and lateral lattice geometry designed to provide robust radiographic imaging unimpeded by reducing overall titanium material and device density.<sup>6</sup>

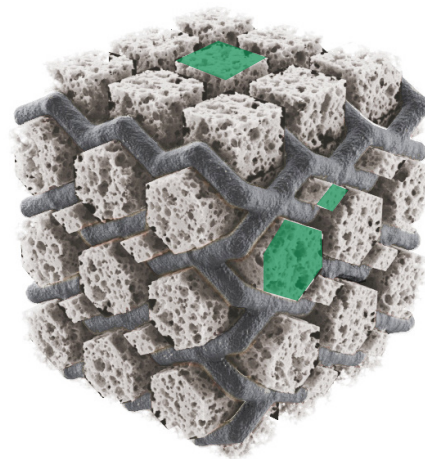

*Image represents potential volume for bony in-growth*

### Studies referenced for the foundational design of NEXXT MATRIXX®

1. Karageorgiou V, Kaplan D. Porosity of 3D biomaterial scaffolds and osteogenesis. *Biomaterials*. 2005;26(27):5474–91.
2. Olivares-Navarrete R, Hyzy SL, Slosar PJ et al. Implant materials generate different peri-implant inflammatory factors: poly-ether-ether-ketone promotes fibrosis and microtextured titanium promotes osteogenic factors. *Spine*. 2015;40(6):399–404.
3. Olivares-Navarrete R, Hyzy SL, Gittens RA, et al. Rough titanium alloys regulate osteoblast production of angiogenic factors. *Spine J*. 2013;13(11):1563–70.
4. Ponader S, von Wilmsowky C, Widenmayer M, et al. In vivo performance of selective electron beam-melted ti-6al-4v structures. *J Biomed Mater Res A* 2010;92A:56–62
5. Li JP, Habibovic P, et al.: Bone ingrowth in porous titanium implants produced by 3D fiber deposition. *Biomaterials* 28:2810, 2007.
6. Data on file at Nexxt Spine, LLC.




# NEXXT MATRIXX®

3D Printed Porous Titanium


## STAND ALONE ALIF

Simple, Integrated, Secure


### PRODUCT FEATURES

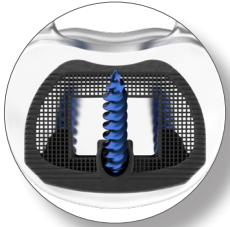


**Anatomically matched profile**  
for appropriate endplate coverage  
and placement on apophyseal rim  
for stability

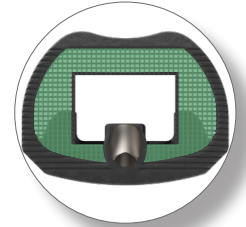

**Ample graft window**  
balanced with lattice landscape  
to create environment for  
bone growth

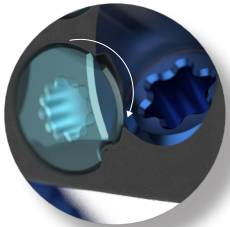


**Integrated one-step turn lock**  
feature to prevent Screw backout

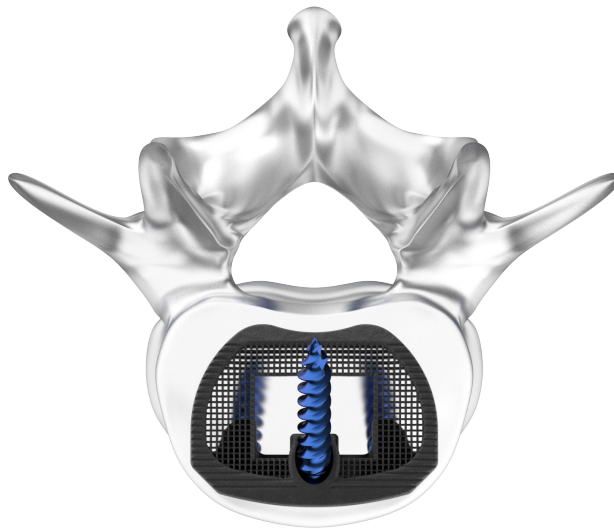

**Self-tapping Screws**  
designed with tip-to-head  
thread pattern for cancellous  
and cortical bone fixation

**Optimized location of Screw Pockets**  
to allow for consistent bone purchase  
for each Screw regardless of Cage height.

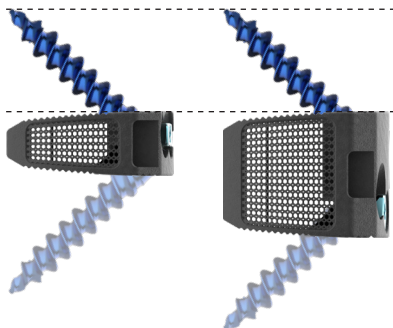

**Multiple Insertion Instrumentation Options**  
to accommodate varying patient anatomies.







# NEXXT MATRIXX®

3D Printed Porous Titanium


## STAND ALONE ALIF

Simple, Integrated, Secure

### CAGE SPECIFICATIONS


#### Footprints

24D x 32W, 27D x 36W, and 30D x 40Wmm

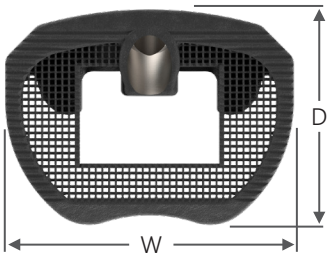


#### Lordoses

8°, 14°, 20°, and 25°\*

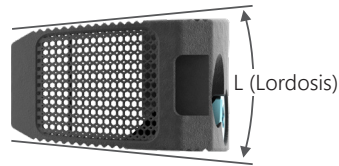

#### Heights


10, 12, 14, 16, 18, and 20mm\*

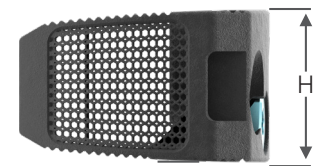

### SCREW SPECIFICATIONS


#### Lengths

20, 25, 30, and 35mm

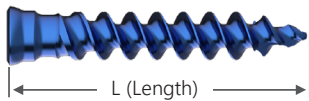

#### Diameters

5.0 & 5.5Ø

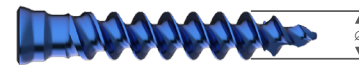

#### Size Offering

Key: Standard ☒ Made to Order ☒

#### Cage

Depth x Width x Height x Lordosis

D	W	H	8°	14°	20°	25°
24	32	10	✓	✓		
		12	✓	✓	✓	
		14	✓	✓	✓	✓
		16	✓	✓	✓	✓
		18	✓	✓	✓	✓
		20	✓	✓	✓	✓
27	36	10	✓	✓		
		12	✓	✓	✓	
		14	✓	✓	✓	✓
		16	✓	✓	✓	✓
		18	✓	✓	✓	✓
		20	✓	✓	✓	✓
30	40	10	✓	✓		
		12	✓	✓		
		14	✓	✓	✓	✓
		16	✓	✓	✓	✓
		18	✓	✓	✓	✓
		20	✓	✓	✓	✓

#### Screws

Diameter x Length

Ø	L	
5.0	20	✓
	25	✓
	30	✓
	35	✓
5.5	20	✓
	25	✓
	30	✓
	35	✓

**NEXXT**  
Spine  
GROW WITH US

**Nexxt Spine, LLC**  
14425 Bergen Blvd, Suite B  
Noblesville, IN 46060  
(317) 436-7801  
Info@NexxtSpine.com  
NexxtSpine.com



For indications, contraindications, warnings, precautions, potential adverse effects and patient counseling information, see the package insert or contact your local representative; visit NexxtSpine.com for additional product information.

All rights reserved. All content herein is protected by copyright, trademarks and other intellectual property rights owned by Nexxt Spine, LLC and must not be redistributed, duplicated or disclosed, in whole or in part without the expressed written consent of Nexxt Spine, LLC. This material is intended for health care professionals, the Nexxt Spine sales force and authorized representatives. Distribution to any other recipient is prohibited.